


Control Box CB4P


CB4P control box is mainly designed for the medical care applications, such as hospital bed and nursing bed. With the optional High Power and Synchronization functions, there are three versions in this series as CB4P, CB4P-SY and CB4P-HP.

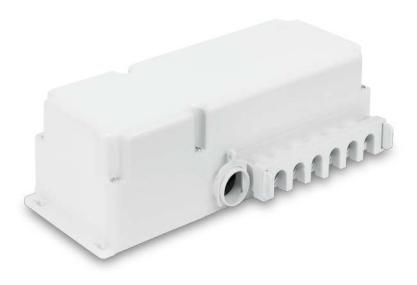


Features and Options

- Main applications: Medical, Home care
- Input: 110V AC 60Hz / 230V AC 50Hz
- Output voltage: 24V DC
- Output rated current: Refer to Compatibility (p.4)
- Power supply: Toroidal transformer 156W / 180W / 216W
- Max. number of actuators: 4 channels
- Max. number of controls: 2 channels
- Duty cycle: 10%, max. 2 min. continuous operation in 20 min.
- Over-current protection function
- IP level: IP54, IPX6
- Electric shock protection level: Class I
- Color: Light grey RAL 7035
- 3-Pin detachable AC power cord: 3000mm
- Rechargeable lead-acid battery (1.3 Ah)

Dimensions

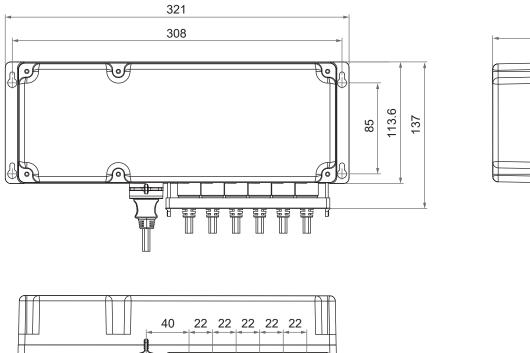


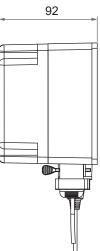


Unit: mm

I/O Introduction


Standard (CB4P)

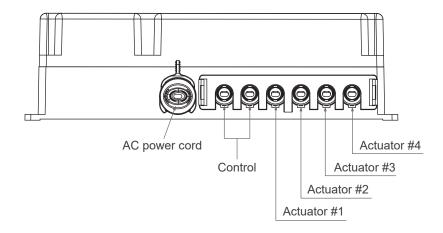

Synchronization version (CB4P-SY)

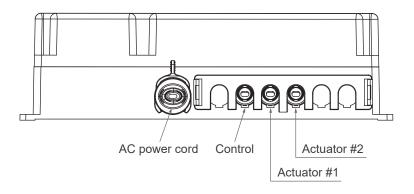

High power version (CB4P-HP)

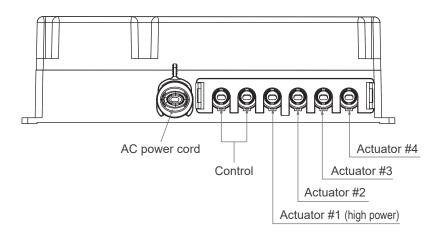

3

Compatibility

Actuators

Version	Actuator port	Compatible models	Compatible requirement
CB4P (Standard)	Max. 4 channels with max. 3A current per channel.	MD66, MD7, MD21, DLC1, DLC2	Without positioning sensor With Moteck H-type DIN plug
CB4P-SY (Synchronization)	2 channels with max. 4.5A current per channel.	MD60, DLC1	Dual Hall effect sensors With Moteck H-type DIN plug
CB4P-HP (High power)	1 channel with max. 9A current.	MD80, MD100	Without positioning sensor With Moteck J2-type phone jack
	Max. 3 channels with max. 3A current per channel.	MD66, MD7, MD21	Without positioning sensor With Moteck H-type DIN plug



The feedback signal phase of dual Hall effect sensors from actuators must be "Type A" defined by Moteck in order to be compatible with CB4P-SY.

Controls

Version	Control port	Compatible models	Compatible requirement
CB4P (Standard)	Max. 2 channels	HC, HF, HT, H3F, H3T, HS12, FS2	Moteck P1 control category With Moteck H-type DIN plug
CB4P-SY (Synchronization)	Max. 1 channel	НС	Moteck S1 control category With Moteck H-type DIN plug
CB4P-HP (High power)	Max. 2 channels	HC, HF, HT, H3F, H3T, HS12, FS2	Moteck P1 control category With Moteck H-type DIN plug

Socket

Control socket for Moteck H-type DIN plug

- CB4P (Standard): Moteck P1 control category
 - 3 drives model / 4 drives model
- 2 drives model







5-pin socket

• CB4P-SY (Synchronization version): Moteck S1 control category




5-pin socket




Actuator socket for Moteck H-type DIN plug

• Without positioning feedback sensors

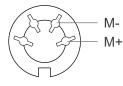


4-pin socket

• Positioning feedback with Hall effect sensor x 1 or x 2

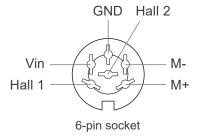

Actuator socket for J2-type Phone jack

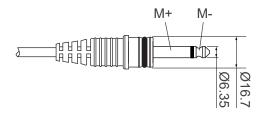

Note: Connect M+ to "Vdc +" & M- to "Vdc -" of DC power to extend the actuator. Switch the polarity of DC input to retract it.



AC Power cord

Plug type (plug in the wall)	Connector type (plug in the control box)	Remarks
US 3-pin plug	IEC 60320 C13 type connector	AC power cord length: 3000mm straight Color: Light grey RAL 7035
CEE 7/7(VDE 3P) plug		
UK 3-pin plug		
CHINA CCC 3-pin plug		

Ordering Key

